SELECIONE O IDIOMA

SELECT THE LANGUAGE

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
Elastic Ligature

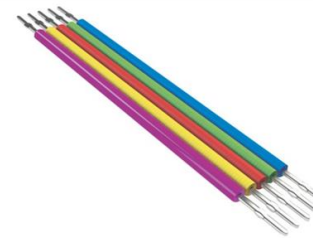They are orthodontic elastomers, used for fixing orthodontic wires or arches to the bracket slots, promoting the transfer of the forces applied by the arches / wires to the teeth.

Use Indication:

- Fixing the arch wire to the bracket.

Handling:

The insertion and removal of the ligatures must be carried out by means of tools to this purpose such as: **Elastic Ligature Applicators** (Ref. 75.01.001) (Ref. 75.01.002) (Ref. 75.01.003) (ref. 75.01.004) or **Tweezer "Plic"** (ref. 75.01.028).

Contraindications:

Morelli supplies the products only to qualified professionals. It is the responsibility of the orthodontist to identify conditions that may be contraindicated for treatment, such as:

- Patients with poor oral hygiene;
- Patients unable to cooperate with treatment;
- Patients with propitious oral environment to dental enamel demineralization;
- Allergy to any of the materials of the orthodontic appliance;
- Diseases or other pre-existing conditions that may hinder orthodontic treatment;
- Root resorption;
- Existing bone resorption;
- Existing Decalcification in dental enamel;
- Use in a Magnetic Resonance environment.

Side effects:

It is the responsibility of the qualified professional to identify any possibility of risk of harm and / or contraindication that may arise during treatment, communicate to the patient the possibility of any unwanted side effect.

During treatment, unwanted side effects may include: damage to oral mucosa and tongue, difficulty in speaking or chewing, tooth discoloration, decalcification, root resorption, periodontal complications, allergic reactions, difficulties in maintaining oral hygiene, discomfort, pain, discomfort and sensibility.

Expiry date:

3 years

Warnings:

Store in a cool and airy place. Avoid exposure to high temperature or intense light, in order to avoid the decrease of the useful life of the product and loss of the mechanical properties.

The product follows the dimensional patterns, initial strength (F_0) and residual after 24h (F_{24}), as established in the norms ABNT NBR ISO 21606 and DIN EN ISO 21606.

Keep the product in its original packaging and always closed in order to avoid drying.

Elastomeric ligatures are inherent to color change in the oral environment: Consider replacing them for shorter periods to maintain the desired aesthetic effect.

This product must not be used in a Magnetic Resonance environment as it may cause interference with the images during the examination. It is up to the orthodontic practitioner to warn the patient about this contraindication and arrange for removal of the device if necessary.

Non reusable product, as its reuse may cause cross infection and loss of mechanical properties due to natural wear.

The product is supplied in unsterile condition, however, has proven biocompatibility, not being necessary cleaning, disinfection or sterilization measures, since the package has not been violated. If the product comes into contact with surfaces or substances other than its intended purpose prior to use, it is recommended to discard the product. The product should not be reused or reprocessed.

Do not use the product if its packaging has been damaged or opened, in order to minimize possible dangers related to the transmission of microorganisms or possible falsification of products.

Product for exclusive application by a qualified professional, the use by person without the necessary technical knowledge may cause unwanted tooth movement, bone loss, loss of the dental element, fenestration, gingival recession or root dehiscence.

It is up to the qualified professional to guide the patient about correct oral hygiene, to avoid the appearance of plaque and tartar, or diseases such as gingival inflammation (gingivitis), periodontitis or even endocarditis.

Considering that the product is used in contact with mucous membranes and body fluids, it is recommended to the professional to apply the sanitary norms applicable for the disposal of medical products according to current sanitary regulations.

ATTENTION:

Morelli is not responsible for the results obtained by applying the product without due precautions or non-observance of the warnings.


Specifications according to DIN EN ISO 21606:

Internal Diameter (D_i); External Diameter (D_o); Cutting thickness in cross section (t); Initial force (F_0); Residual force after 24 hours (F_{24})

Elastic Ligature in Canes ($D_i - \varnothing 1,4^{+0,05/-0,05}$); ($D_o - \varnothing 3,0^{+0,1/-0,1}$); ($t - 0,07^{+0,05/-0,05}$); ($F_0 - 350gf \pm 150gf / 3,4N \pm 1,5N$); ($F_{24} - \rightarrow = 50\%$)

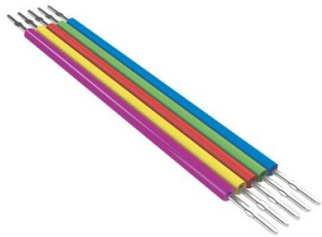Elastic Ligature in Modules ($D_i - \varnothing 1,4^{+0,1/-0,1}$); ($D_o - \varnothing 3,1^{+0,05/-0,05}$); ($t - 0,85^{+0,05/-0,05}$); ($F_0 - 350gf \pm 102gf / 3,4N \pm 1N$); ($F_{24} - \rightarrow = 45\%$)

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